



Specialized Technology Resources
 85 John Road
 Canton, MA 02021
 Phone: 781-821-2200
 Fax: 781-821-9266
 www.STRQuality.com

Technically Advanced Quality Assurance TAQA® Evaluation Report

Mr. Mahmoud M. Sales
SCI International, Inc.
 5902 Enterprise Court
 Fredrick, MD 21703

Folder #: 1013383-1 & 1015170
 Date: 6/25/2010 & 7/23/2010
 Received Date: 6/8/10 & 7/1/10
 Page: 1 of 3

cc: Mr. Stuart Neuman

Product:	SCI Brand Perfect Choice Pregnancy Midstream Test*	UPC Code:
Size:	Sample	Reference #: Prepaid American Express
Supplier:	SCI International, Inc.	Inventory ID #: 239461 & 242896*
Code:	Lot 06119M11 Exp 12/2011 & Lot 12159M11 Exp 07/2012*	Lab Sample #: 1013383001 & 100517002*
Purpose:	TAQA EVALUATION vs. FIRST RESPONSE	

National Reference Brand: First Response Pregnancy Test Early Result (6 days before a missed period);
 Lot Code: BU0068PA Exp 02-2012 - provided by SCI International

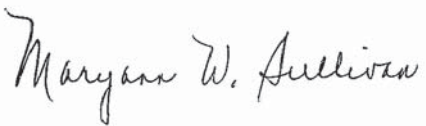
Overall Rating: 10.0 Sample submitted meets TAQA quality criteria for parameters examined (a).

Comments: (a) Differs in construction from First Response but result/performance is comparable.

Note: * SCI Brands includes: New Choice, Perfect Choice, Baby Test, Baby Check, UCheck

Overall Quality Rating *

- 9.0 - 10.0** Sample tested is Comparable to the National Brand; meets TAQA quality criteria
- 7.5 - < 9.0** Sample tested is Marginally Comparable to the National Brand; does not meet all TAQA quality criteria
- < 7.5** Sample tested is Not Comparable to the National Brand; does not meet TAQA quality criteria

 Digitally signed by Maryann
 Sullivan
 Date: 2010.07.23 13:04:02 -04'00'

Maryann W. Sullivan, Client Services Manager

Determination		First Response Pregnancy Test Early Result (6 days before a missed period)	SCI Brand Perfect Choice Pregnancy Midstream Test	Comments
Performance - Bench Top Used Sample Inv.#242896	Positive Results Used HCG Hormone Positive Control 20miu/mL Lot 060710P Exp 6/20/10 (supplied by SCI)	Using a 5 second dip "urine" entered result window at 56 seconds; first reaction at 75 seconds; complete results at 80 seconds; constant for 3 minutes	Using a 5 second dip "urine" entered the window at 12 seconds; first sign of a reaction at 28 seconds; complete positive result at 38 seconds; constant for 10 minutes	
	Negative Results Used Quantimetrix The Dipper Lot: 44501 Exp 2011-01	Using a 5 second dip "urine" entered result window at 51seconds; first reaction at 75 seconds; complete results at 85 seconds; constant for 3 minutes	Using a 5 second dip "urine" entered the window at 12 seconds; first sign of a reaction at 30seconds; complete positive result at 40 seconds; constant for 10 minutes	
Physical	DIMENSIONS, inches Overall Length Width at Cap Width at Handle	5-5/8 13/16 3/4	5-5/8 5/8 3/8 (at tip)	
	ABSORBENT TIP Length Width Thickness	1-1/16 7/16 1/8	1-3/4 3/8 1/8	
	HANDLE Length Width Thickness	4-3/8 13/16 5/16	3-1/2 5/8 - 3/8 (tapered) 5/16	
	CAP Length Width Thickness	1-5/8 1-3/16 5/16	2-1/8 (tapered) 5/8 - 3/8 (tapered) 5/16	

		First Response Pregnancy Test Early Result (6 days before a missed period)	SCI Brand Perfect Choice Pregnancy Midstream Test	Comments
Physical	Construction	Opaque/white plastic contoured handle with clear/pink plastic tapered cap; white rectangular-shaped absorbent tip at one end. At mid section, one oval result window covered with plastic with printed results and symbols to one side of it; thumb grips at end of handle; First Response and symbol embossed on lower handle.	Opaque/white plastic tapered stick with purple tapered plastic cover; stick has an extruded (plastic) sampling collector "tray" containing an absorbent strip at one end, a rectangular result view window at mid section and a depressed finder hold at the other end. (a) (samples submitted)	(a) Differs in construction from First Response but results/performance is comparable
	Product Workmanship	Acceptable	Acceptable	
	Packaging	Individually wrapped in printed heat-sealed foil pouch; boxed with clear plastic overwrap. Instruction pamphlet with illustrations inside.	Individually wrapped sticks in heat-sealed foil packets (use instructions provided)	